


UniBloc Analytical Balances

AUW-D/AUW/AUX/AUY Series






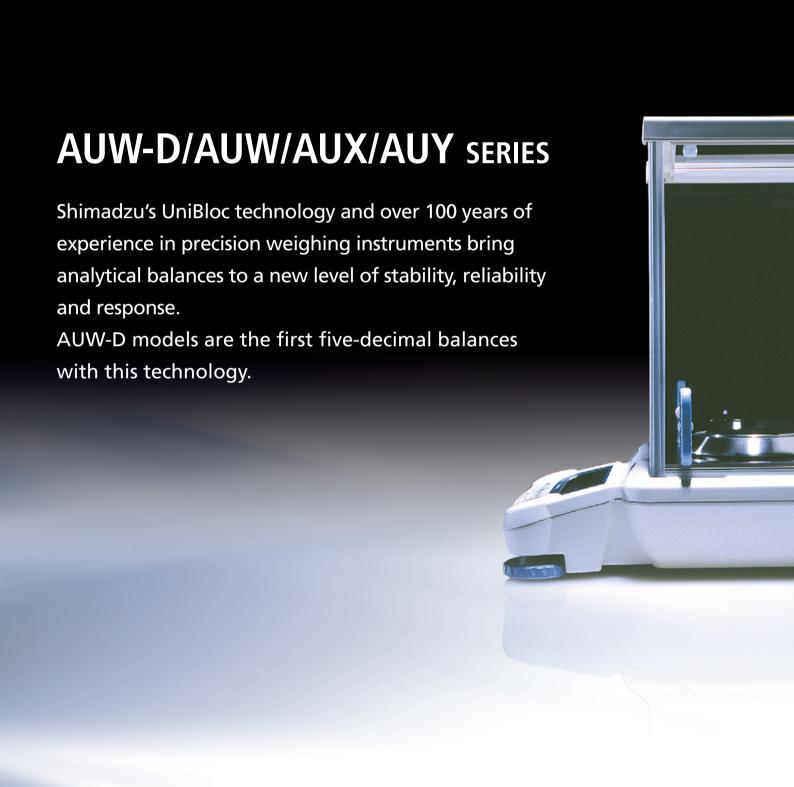





Shimadzu introduced one piece force cell technology for precision balances in 1989. The Today's UniBloc is created by high-precision electric discharge wire processing applied to a block of aluminum alloy, and replaces the conventional electro-magnetic balance sensor assembly. UniBloc's compact, uniform structure ensures stable temperature characteristics, excellent response time and stable corner-load performance. The UniBloc design permits a consistancy of production that assures reliability and a long operational life. AUW-D dual range analytical balances are the first five-decimal balances with the advantages of UniBloc one-piece force cell technology.

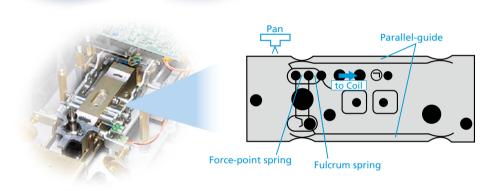



Excellent Weighing Performance

Fast Response

Compact UniBloc mechanism and digital processing technology produce fast response and stability at the same time.

Environmental adjustment

Microprocessor digital control can be set to automatically provide the most suitable data processing for the installation environment and weighing application.

Measurement Administration



★ GLP/GMP/ISO9000 Conformance (AUW-D/AUW/AUX models)

If an optional printer is connected, data can be printed out with the date and time. Sample ID and counting up numbers can also be attached to the data with EP-110 printer. Calibration report can be automatically output, ensuring the measuring control and traceability required by GLP/GMP/ISO9000.



DATE 2018-02-31 TIME 14.45.37 ID: 780315 ID: 780315 No.0010203001 0.074029 DATE 2018-02-31 TIME 14.46.11 ID: 780315 No.0010203002 0.04959s DATE 2018-02-31 TIME 14.46.39

Automatic Span Calibration

Span calibration is essential for accurate measurement in precision weighing to adjust for the effects of even small changes in ambient temperature. AUW-D/AUW/AUX models have a built-in motor driven calibration weight to automate this necessary operation. The balance can be set to take care of the calibration itself, leaving the operator to concentrate on measurement work.



Fully-automatic calibration by temperature detection: PSC

The balance detects variations in the ambient temperature that influence accuracy and automatically performs calibration to compensate it. (AUW-D/AUW/AUX models) This great feature has been provided by Shimadzu since as early as 1985.



Fully-automatic calibration at user pre-set times: Clock-CAL

The balance automatically performs calibration at selected times, up to three times a day (e.g., before starting work, during lunch-break, or after work). (AUW-D/AUW models)





Span calibration at any time: **Touch-key calibration**

Automated calibration can be started by pressing keys. (AUW-D/AUW/AUX models)

Also, your external calibration weights can be used for span calibration. (All models)



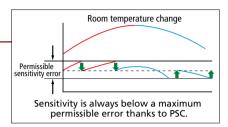

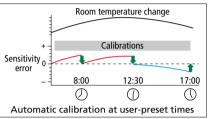

Application Support



Efficient filling or dosing work

"Pouring mode" brings very fast response with even a very small amount of sample added. Suitable for filling or dosing works.



Formulation mode

Convenient for making many measurements of minute samples and seeking the total mass.



Specific Gravity Measurement

Installing the optional SMK-401 specific gravity kit transforms the balance into a dedicated instrument for measuring specific gravity or density. Specific gravity measurement software is already installed in the Shimadzu balance.



Below-weigh Hook for Hanging Measurement

Measurement which cannot be performed within the weighing chamber is possible using the hook provided as a standard component.



Piece Counting and Unit Conversion

In addition to piece counting, the balance can also perform weight measurement as percentages and in a variety of mass units, such as carat.



Interval Timer

Data can be automatically output at time intervals set in the range from 1 second to 99 minutes 59 seconds. (AUW-D/AUW/AUX models)





* with optional AKB-301 Application Keyboard



User-friendly Features



User-friendly Weighing Chamber

The large dimensions of the weighing chamber allow easy measurement of even tall flasks. Equalization in temperature before measurement is essential in precision weighing. This is easy with the extra space allowing samples to be kept inside the chamber. Weighing work is made easy by the smooth door movement. The doors can be detached to allow the chamber to be cleaned with ease and it is also possible to remove, clean or replace the door rail.







Easy Level Adjustment

Level adjustment, particularly important in installing an analytical balance, can be performed with ease using an easy-to-view level gauge at the front of the balance and large level screws.





Metal Housing

The aluminum die-cast housing offers lasting high quality exterior, protection of the core mechanism, and ease of cleaning.





Backlit Display

Allows use under poor lighting conditions (AUW models)



Comfortable Key Operation

The embossed key panel sheet provides clear clicking response when operated. If you choose, the key operations are confirmed with a gentle beeping sound, too.



💂 Analog Bar Graph Display

The analog display allows the operator to see at a glance how much more can be loaded before reaching the weighing capacity.




Keeps dirt away from key panel and display.

Balance Keys

Data directly jumps to your Windows® application

Shimadzu's unique offer: the weighed result on the balance's display is "typed" at the cursor position or activated cell in any application on Windows. This advantage requires no software installation. A cable is all that you need.

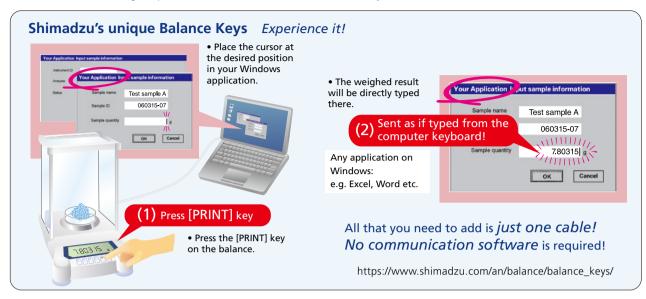

Expanding Possibilities with a Wide Range of Accessories (Optional)

Name	Specification
2-way Ionizer (Static Remover) STABLO-AP	Secure static removal; Air blower switched ON/OFF.
Electronic Printer EP-100	Impact-dot print
Electronic Printer EP-110	Impact-dot print; Statistical calculation; Sample No. can be attached; Date/time can be attached. (except for AUY)
Specific Gravity Measurement Kit SMK-401	This kit makes it possible to determine specific gravity or density of solids with the highly precise immersion method. Samples with weights up to the weight capacity of the balance can be measured.
RS-232C Cable	Required to connect to a PC.
Application Keyboard AKB-301	Enables unit weight settings for piece counting, unit weight display, recalculation of piece with new unit weight, taring and output of currently displayed value. Other numerical settings can also be performed more efficiently.
Foot switch FSB-102PK	For hands-free print command (key connector)
Foot switch FSB-102TK	For hands-free tare command (key connector)









AKB-301

Consumables and Replacement

Name	Specification					
In-use Protective Cover (standard accessory)	Cover for display (1 cover provided as standard accessory; 5 covers in a set)					
Printer paper for EP-100 / EP-110	(20 rolls)					
Ink ribbon cartridge for EP-100 / EP-110	Standard type, black ink (5 pieces)					
Ink ribbon cartridge for EP-100 / EP-110	Long life type, black ink (5 pieces)					









Specifications

	Series	AUW-D (dual range semi-micro)		AUW			AUX			AUY		
	Model	AUW220D	AUW120D	AUW320	AUW220	AUW120	AUX320	AUX220	AUX120	AUY220	AUY120	
Capacity		220 g / 82 g	120 g / 42 g	320 g	220 g	120 g	320 g	220 g	120 g	220 g	120 g	
Minimum display		0.1 mg / 0.01 mg	0.1 mg / 0.01 mg	0.1 mg	0.1 mg	0.1 mg	0.1 mg	0.1 mg	0.1 mg	0.1 mg	0.1 mg	
Repeatability (Standard deviation, σ)		≤0.1 mg (large range) ≤0.05 mg (small range)	≤0.1 mg (large range) ≤0.02 mg (small range)	≤0.15 mg	≤0.1 mg	≤0.1 mg	≤0.15 mg	≤0.1 mg	≤0.1 mg	≤0.1 mg	≤0.1 mg	
Linearity		±0.2 mg (large range) ±0.1 mg (small range)	±0.2 mg (large range) ±0.03 mg (small range)	±0.3 mg	±0.2 mg	±0.2 mg	±0.3 mg	±0.3 mg	±0.2 mg	±0.2 mg	±0.2 mg	
Response time (typical)		3 s (large range) 15 s (small range)	3 s (large range) 12 s (small range)									
Operating ambient temperature		5 to 40 °C										
Temperature coefficient of sensitivity (10 to 30°C)		±2 ppm/°C (When PSC is Off)							±2 ppm/°C			
Sensitivity stability against temperature change (When PSC is On, 10 to 30°C)		±2 ppm								-		
Pan size		80 mm dia. approx.										
Body dimensions		W220 mm × D330 mm × H310 mm approx.										
Weight		7 kg approx.										
Power consumption		7 VA approx.										
	Backlighted display				•	•						
	Built-in calibration weight	•	•	•	•	•	•	•	•			
	PSC	•	•	•	•	•	•	•	•			
Functions	Clock-CAL	•	•	•	•	•						
	Built-in clock	•	•		•	•	•	•	•			
	GLP/GMP/ISO calibration report	•	•	•	•	•	•	•	•			
and	Balance Keys	•		•	•	•	•	•	•	•	•	
features	Formulation mode	•	•	•	•	•	•	•	•	•	•	
	Interval timer output	•	•	•	•	•	•	•	•			
	RS-232C I/F	•	•	•	•	•	•	•	•	•	•	
	Specific gravity measurement software, piece counting, % display	•	•	•	•	•	•	•	•	•	•	
	Analog display	•	•	•	•	•	•	•	•	•	•	

^{*} Windows, Excel and Word are registered trademarks of Microsoft Corporation.



Shimadzu Corporation www.shimadzu.com/an/

For Research Use Only. Not for use in diagnostic procedures.
This publication may contain references to products that are not available in your country. Please contact us to check the availability of these

Inis publication may contain references to products that are not available in your country. Please contact us to check the availability of these products in your country.

Company names, products/service names and logos used in this publication are trademarks and trade names of Shimadzu Corporation, its subsidiaries or its affiliates, whether or not they are used with trademark symbol "TM" or "®".

Third-party trademarks and trade names may be used in this publication to refer to either the entities or their products/services, whether or not they are used with trademark symbol "TM" or "®".

Shimadzu disclaims any proprietary interest in trademarks and trade names other than its own.

The contents of this publication are provided to you "as is" without warranty of any kind, and are subject to change without notice. Shimadzu does not assume any responsibility or liability for any damage, whether direct or indirect, relating to the use of this publication.